


Made from 100% food-safe silicone

All products that come into contact with your body and precious breast milk should meet the highest quality standards. The KindestCup has therefore been tested according to FDA and LFGB standards as a product with food contact and found to be safe. Thanks to the large opening, the KindestCup is easy to clean by hand and consists of just one single piece.

Universally suitable for every breast size

Due to the large opening of the KindestCup, the suction effect is not only limited to the nipple, but also acts on the surrounding skin. Therefore, neither the size of the nipple nor the size of the breast plays a role in the suction effect of the KindestCup.

KindestCup

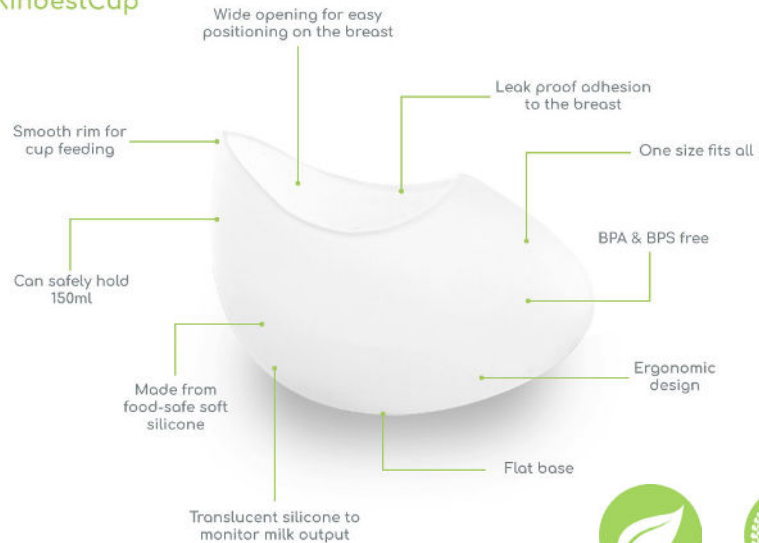

BPA free



WHO code compliant

ardo
SWITZERLAND



KindestCup Milk Collection Cup

At Ardo we understand what it is like to be a new mum and to feel overwhelmed by having a newborn baby and needing to feed.

This leaflet illustrates how to use the KindestCup



SCAN ME


View our range of expressing and information videos here

ardo
SWITZERLAND

Manual Expressing

The KindestCup is designed specifically for gently obtaining breast milk with active manual expression (compressing the breast with your hand)


- The large opening makes it easier to position the nipple and ensures that not a drop of breast milk is lost during manual expression.
- The soft, flexible silicone fits comfortably in the hand, is easy to grip and enables convenient one-handed use.
- The sturdy base provides firm support - ensuring that valuable breast milk is collected reliably and spillage prevented.

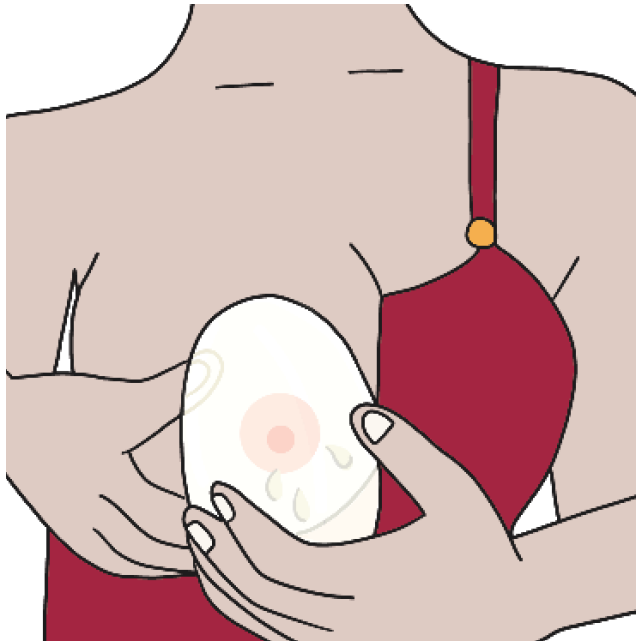

Collection Cup

The KindestCup is a collection cup for breast milk which flows on its own thanks to the let-down reflex

- Here, the baby breast feeds or the breast milk is actively expressed on one breast, while the KindestCup stays on its own on the other breast thanks to a gentle vacuum and collects flowing milk which would otherwise be lost in a nursing pad.
- Even lying on your side while breastfeeding at night, the KindestCup stays leak-proof and gives the mother a feeling of security.
- The KindestCup does not actively encourage breast milk expression; the vacuum simply provides secure grip on the breast. Unlike conventional silicone breast pumps, the KindestCup therefore does not promote over production, but supports on-demand milk production in the natural breastfeeding rhythm.

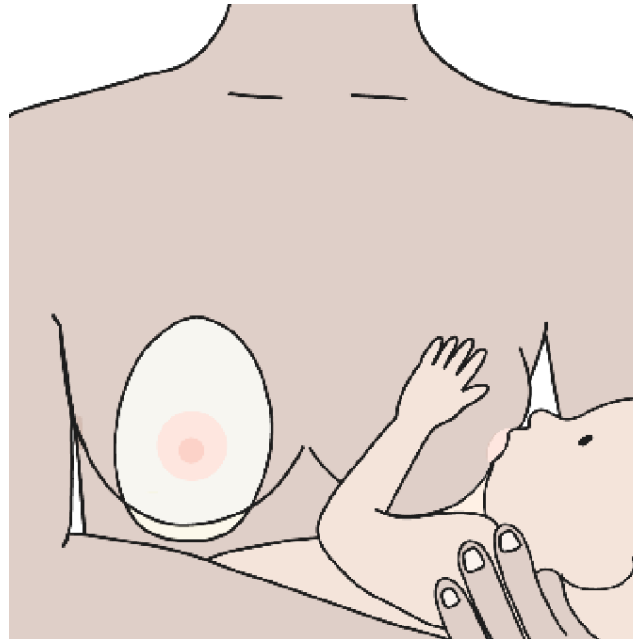

Feeding Cup

The KindestCup is also a feeding cup

- Here, the KindestCup is held in one hand and the edge moved to the baby's mouth. The baby starts to drink the milk on its own once it comes into contact with the KindestCup.
- If the milk is collected directly with the KindestCup and, in turn, fed to the baby, no milk is lost by having to transfer the milk into another container. That is particularly useful when obtaining and feeding colostrum to the baby.

