



SUCTION SENATOR & MASTER

- ✓ reliable
- ✓ powerful
- ✓ very quiet



SWISS QUALITY 

 ardo
SWITZERLAND

VACUUM PUMPS

RELIABLE, POWERFUL, AND QUIET



Both the Master and Senator vacuum pumps impress users thanks to their uncompromising suction power, rapid vacuum generation, and many other benefits:

- ✓ Large vacuum meter simplifies setting and reading the vacuum.
- ✓ High-quality drive unit with glass cylinder and graphite pistons guarantee the highest level of precision and reliability with extremely quiet operation.
- ✓ Hydrophobic bacteria filter for optimal safety.
- ✓ Foot-operated vacuum regulator and foot switch available for hands-free operation.
- ✓ Easy-to-clean surfaces for simple hygiene.

The user-friendly pumps with their modern design provide users with optimal support for a wide range of procedures and are silent helpers in the background.

MASTER

Suction power 50 l/min,
vacuum up to 95 kPa, 713 mmHg



SENATOR

Suction power 30 l/min,
vacuum up to 95 kPa, 713 mmHg



incl. fitting rails

Master & Senator
also available with
foot switch

for hands-free switching on/off



The Master is the high performance version among the Ardo vacuum pumps with uncompromising suction power of 50 l/min.

Approved areas of application:


- ✓ Surgery
- ✓ Liposuction
- ✓ Dental and maxillary surgery
- ✓ Veterinary medicine

The Senator is the right choice for all applications with a suction power of up to 30 l/min.

Approved areas of application:

- ✓ Endoscopy
- ✓ Gynecology & obstetrics
- ✓ Ear, nose, throat (ENT)
- ✓ General suction of secretions

Ardo pumps are very quiet and extremely low maintenance.

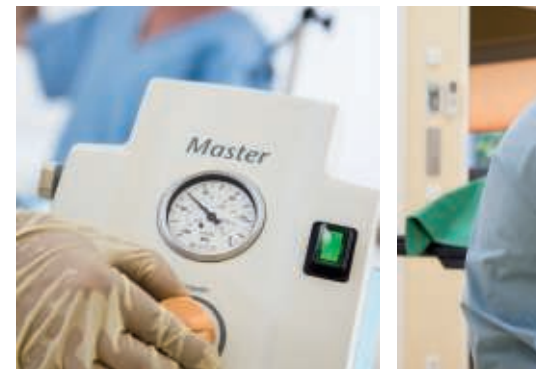

OUR BENEFITS

WHY CHOOSE AN ARDO VACUUM PUMP?



3 guard system

Optimal safety

The "3guard" safety system provides optimal protection and the highest level of safety in all situations. The hydrophobic bacteria filter has a double function. It protects the aggregate in case of possible overflow of liquid and foam while also providing highly effective protection against bacteria. The float integrated into the lid of the suction jar effectively prevents a full suction jar from overflowing.



LowFriction

High-quality piston system

The Pyrex glass cylinder with the corresponding graphite pistons forms the heart of the unique "LowFriction" cylinder/piston system. The self-lubricating, maintenance-free system promises almost silent operation for optimal concentration during application.

The combination of high-quality components and Swiss precision in production ensure excellent functioning. The rapid vacuum generation and high suction power ensure optimal performance for any type of application.



Reusable accessories

Choose from our reusable accessories according to your needs. Reusable suction jars are available in sizes 1, 2 and 5 litres, made of high-quality polysulfone, autoclavable, and have easy-to-read graduations. The lid fits all jar sizes and the holder for various fitting rails is already integrated in the lid.



Individual application sets

Assemble your own equipment or select from various pre-configured sets. With the many options available from Ardo, you will always be perfectly equipped for any occasion.

- ✓ Durable and robust trolley
- ✓ Foot-operated vacuum regulator—with a matching holder for the trolley if required
- ✓ Change-over valve for parallel switching of two suction jars
- ✓ Instrument tray suitable for fitting rails
- ✓ 3 suction jar options



*In delicate situations,
you can rely completely
on Master and Senator
vacuum pumps.*



PRODUCT RANGE

MASTER & SENATOR VACUUM PUMPS AND ACCESSORIES



Master & Senator

Item no.	Description	Voltage	Hertz
30.00.31	Master Vacuum Pump 50 ltr/min	230V	50/60Hz
30.00.32	Senator Vacuum Pump 30 ltr/min	230V	50/60Hz
50.00.170	Stand, Senator/Master		

Accessories for Vacuum Extraction

Item no.	Description
50.00.167	1Ltr suction jar, polysulphone, graduated
50.00.168	2Ltr suction jar, polysulphone, graduated
50.00.169	5Ltr suction jar, polysulphone, graduated
50.00.08	Lid with handle and mechanical overflow safety device, for 8mm suction tubing
50.00.07	Foot operated vacuum regulator
50.00.81	Silicone Vacuum-Extractor Cup 50mm with silicone tubing and angled coupling connector
50.00.80	Silicone Vacuum-Extractor Cup 60mm with silicone tubing and angled coupling connector
50.00.25	Bird Cup Anterior Vacuum-Extractor 40mm
50.00.26	Bird Cup Anterior Vacuum-Extractor 50mm
50.00.27	Bird Cup Anterior Vacuum-Extractor 60mm
50.00.24	Traction Handle for Bird Vacuum-Extractor Cup
50.00.214	Bird Cup Posterior Vacuum-Extractor 50mm with traction handle



Contact us for more information:

Ardo medical Ltd
Unit 1,
Belvedere Trading Estate,
Taunton, TA1 1BH

T+44(0)1823336362
info@ardomedical.co.uk
www.ardomedical.co.uk

SWISS QUALITY

99.04.312 / V01 / 06.2019